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EN

Thank you for purchasing the AOT electric massage roller. We are confident that you will be satisfied with the product. Read the user manual before use and keep it for future reference.

1. Specifications

Product name	AOT electric massage roller
Power	40 W
Dimensions	52,5 x 14,5 x 8,6 cm
Battery	4000 mAh, 7.4 V, USB-C charging
Rotational speed	280-590 RPM, 6 speed settings
Red light wavelength	630 nm

2. Warnings and safe use

This product is an electric massage roller designed for adult self-care, for recovery after exercise, warming up muscles and other muscle care. The product is designed for home use and is not intended for professional use. Read the user manual before you start using the product.


SAFETY RISK: Make sure that hair, jewellery or clothing cannot get caught in the rotating parts of the product.

- Check with a physiotherapist or other healthcare professional that using the product is safe for you.
- Use of the product is generally not recommended if:
 - you are pregnant or recovering from surgery
 - you have a joint prosthesis, pacemaker, internally implanted metal support structures or similar
 - you have a herniated disc, migraine or retinal damage
 - you have been diagnosed with diabetes or epilepsy
- This product is not intended to treat, diagnose or cure any disease.
- This product is not a toy. Keep the product out of the reach of children.
- Do not allow children or persons with reduced physical, functional or decision-making ability to use the product without supervision.
- The product must not be used on the face, head, neck or genital area.
- Clean the product with a dry or slightly damp cloth. Never wet the product.

Safe use of a battery-powered electrical product

- Read these instructions carefully before using the product for the first time and keep them for possible future use.
- Use the product only in accordance with the instructions provided in this user manual.
- Use only the charger supplied with the product or an approved charger.
- The battery must not be disassembled, damaged or exposed to heat.
- The product must not be left charging unattended. Do not continue charging if the product becomes hot.
- If the power cord or charger is damaged, do not use or charge the product.
- Do not expose the product or battery to heat, open flame or liquids during use, storage or charging.
- Do not attempt to disassemble or replace the battery yourself.
- Incorrect handling of the battery may cause explosion, fire or personal injury.
- This product is intended for household use only.
- Do not leave the product in direct sunlight or in hot spaces.
- Charge the product on a dry, level surface. Make sure that the temperature of the charging area is 0–40 °C.
- Unplug the charger from the wall socket when charging is complete.
- Check the battery regularly for swelling or damage.
- Store the product in a dry and cool place (10–30 °C).
- Before long-term storage, the battery should be fully charged. If the product will not be used for a very long time, the battery should be charged regularly to at least 50%.

3. Product parts

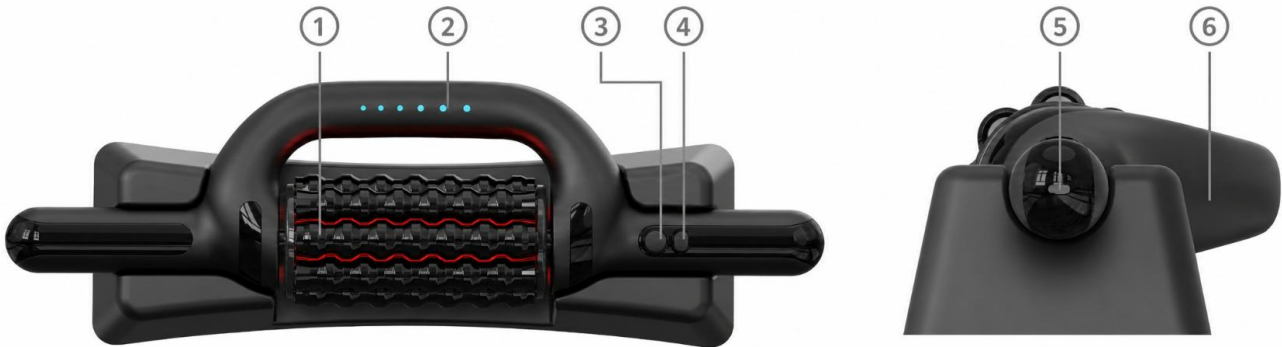1. Massage rollers
2. Indicator lights
3. Power button
4. Red light switch
5. Charging port
6. Handles

4. Using the product

Functions

1. Press the power button for 2 seconds to start the product.
2. Change the power/rotational speed by briefly pressing the button. The product has six (6) speed levels. The indicator lights on the handle show the selected speed (1–6).
3. Change the direction of rotation by double-clicking the power button.
4. Activate the red light using the separate switch next to the power button.

Areas of use

Always make sure before starting the product that hair, jewellery or clothing cannot get caught in the rotating parts of the product.

BACK

CALVES

IT BAND & GLUTES

SHOULDERS

THIGH & SQUADS

HAMSTRINGS

FEET

ARMS

5. Warranty

This product has a 1-year warranty granted by the importer, valid from the date of purchase.

The warranty covers manufacturing and material defects. It does not cover normal wear and tear, or defects caused by user error, damage, incorrect installation or use contrary to the user manual. The warranty will also become void if the product has not been used, maintained or stored with appropriate care. The warranty is valid for household use only.

Keep the receipt or other proof of purchase for warranty handling. We recommend, for example, attaching the receipt or recording any order number or similar identifier in this user manual.

6. Recycling

Dispose of the packaging, packaging materials, product and battery properly.

Do not dispose of the battery or product with household waste. Take it to an appropriate recycling point.

The product must be sorted as electronic waste. Recycling points can be found, for example, in shops that sell electronics and at waste stations.

The product contains a fixed built-in battery, which must only be removed at an appropriate recycling facility. Dispose of the product as electronic waste.